


# 20.01. Vaccine Vial Monitoring (VVM)

- There used to be a comprehensive description of VVM on the WHO website, but, alas, no longer  . We can resort to the [Wikipedia article on VVM](#)
- As of August 2023, WHO has published a [VVM infographic](#). In case that is also removed from the WHO web site, we have stored it here:

VVM infographic


## Configuration

There are a few configuration steps needed to use VVM (Vaccine Vial Monitoring) in your mSupply system.

### Set up vaccine vial monitor (VVM) statuses

To configure VVM status options, you need to have permission in the system to edit vaccine vial monitoring statuses. See the [23.15. Managing users](#) page for details on how to set user permissions. You need to turn on the *View and edit vaccine vial monitor status* permission on the *Permissions (3)* tab.

You can access the vaccine vial monitor (VVM) status configuration window by choosing *Item > Vaccines > Show vaccine vial monitor status...* from the menus. This will open the VVM status window:

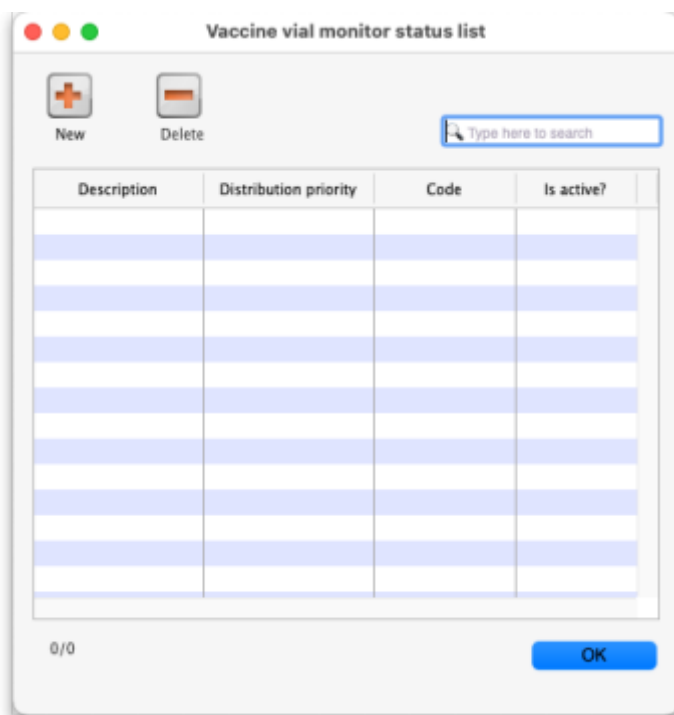

The window will show you a list of all the current vaccine vial monitor statuses that have been setup.


To edit one, simply double-click on it in the list. This will open the same window as when adding a new status (see below) but populated with the status's current settings, ready for you to edit them.

The columns in the table contain the following:

- **Description:** The name of the status. This is what will be displayed for this status throughout mSupply
- **Distribution priority:** The order in which is will be displayed and selected by mSupply for distribution. Lowest number is distributed first. This priority overrides expiry date for vaccines.
- **Code:** The code of the status. A short name by which is can be remembered, not visible elsewhere in mSupply.
- **Is active?:** If checked the status is active. If a status is active it can be assigned to new vaccine stock. Inactive statuses cannot be assigned to new stock but they will be visible where they have been assigned to existing stock or other existing records.

## Adding a VVM status

Click the **New** button. This window opens:

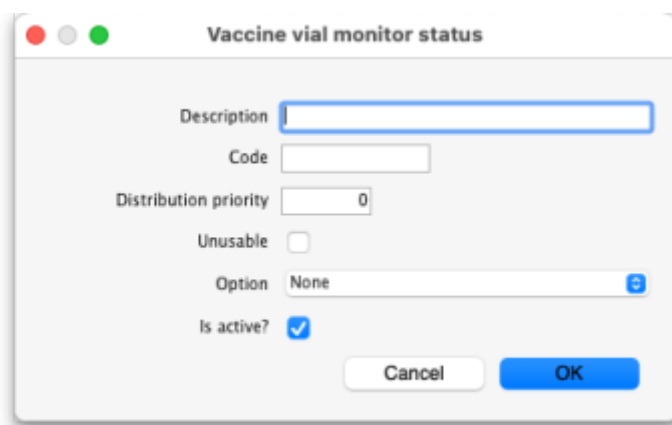

- **Description:** the name of the status and how it will be referred to throughout mSupply.
- **Code:** a short code used to refer to the status.
- **Level:** the numerical level the status belongs to.
- **Is active?:** checked by default. If this is checked then the status will be available for selecting in mSupply mobile.

Click on the **OK** button to save your changes or the **Cancel** button to cancel them.

## Priority of VVM Status

When issuing stock, this will be prioritised based on Level 1 to 4 on mSupply. Level 1 will be prioritised first, then level 2 and so forth. You have the flexibility to customize your VVM stages to determine which should be prioritized first.

In this example, we can see that Stage 2 has first level priority for being issued:




## Deleting a VVM status

You can only delete a status that you have just created. Once this list of VVM statuses has been closed, you cannot delete the status any more. If there is one that you wish to delete but can't then edit it and uncheck its **Is active** checkbox (as described above). That will effectively delete it because it will no longer be selectable in any mSupply mobile functionality (but it will mean the status is still available for historical records).

If the status you want to delete is newly created:

1. Select the status you want to delete by clicking on it in the list
2. Click on the **Delete** button
3. The window below will appear asking you to confirm the deletion. Click **Yes** to confirm.

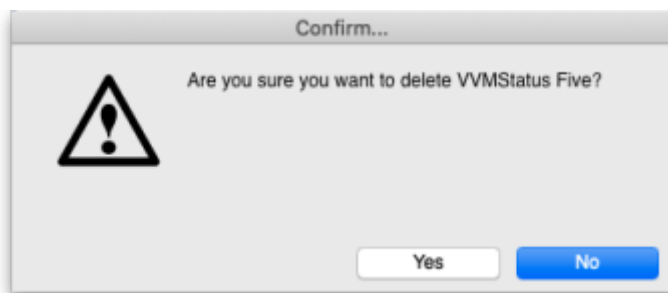

## Tell mSupply which items are vaccines

Tell mSupply which items are vaccines by checking the **Is a vaccine** checkbox and entering the number of doses per vial in the **Doses** field on the *General* tab of the item's details window. See the [4.01.01. Items - General tab](#) page for details.

## Turn on preferences

Finally, you need to turn on some preferences so that VVM status is exposed in the mSupply interface:

- The *Able to specify VVM status when receiving items* store preference to be able to enter the VVM status for a vaccine when it is received. See the [26.07. Virtual stores](#) page for details.
- The *Sort available batches by VVM status rather than Expiry* store preference if you want the VVM status of the stock to be included when deciding what order to offer stock for selection. See the [26.07. Virtual stores](#) page for details.

## Using VVM status

### Enter the VVM status of vaccines when they are received

When the *Able to specify VVM status when receiving items* store preference is turned on (see above) you can enter the VVM status of vaccines when they are received into a store on a supplier invoice. See the [7.01. Supplier invoices](#) page for details.

### Update the VVM status of vaccines when you check them

When you check the status of the VVM labels for a batch, you can enter the details of the current VVM status into the stock line details window.

You can get to this window two ways:

1. On the *Item* tab of the Navigator, click on the *Stock* icon then double-click on the line of stock you are interested in.
2. Search for items, double click on an item to open its details window, click on the *Stock* tab to show the item's current stock and then double-click on the stock line you are interested in.

Whichever way you choose, this window will open:

The screenshot shows a window titled "Item stock line details...". At the top, it displays "Item code bcgvac10" and "Item name BCG Vaccine". Below this are tabs for "General", "Ledger", "Status", "Reference documents", and "Log". The "General" tab is active. The window contains several input fields and labels:

- Quantity according to stock:** Total Quantity 2000, Available 2000
- Cost price:** 2.3
- Sell price:** 2.3
- Expiry date:** 31/08/2026
- Batch:** ABCD1234
- Supplier:** ACME Pharma
- Donor:** None (with a "Change donor" button)
- Note:** A large empty text area.
- Custom stock fields:** A large empty text area.
- Other fields:** Hold (checkbox), Location (None), Pack size (10), Vial, Volume per pack (0 m3), Weight per pack (0 kg), Quan per inner (1), Inners per Outer (1), Initial quantity (0).

At the bottom right, there are "Cancel" and "OK" buttons.

Click on the Status tab:



been entered.

If you enter a VVM status that means the stock is unusable, we recommend that you also check the **Hold** checkbox on the *General* tab (shown above) so that the stock is not issued.

## Check the VVM status of stock when issuing

When the *Able to specify VVM status when receiving items* store preference is turned on (see the [Turn on preferences](#) section above) you will be able to see the current VVM status of vaccine stock and issue in packs or doses when issuing a vaccine item. See the [8.01. Issuing goods to a customer \(customer invoices\)](#) page for details.

Previous: [20. Cold Chain Equipment](#) | | Next: [20.02. Configure VVM for a Mobile store](#)

From:  
<https://docs.msupply.org.nz/docs/> - **mSupply documentation wiki**

Permanent link:  
[https://docs.msupply.org.nz/docs/cold\\_chain\\_equipment:configure?rev=1768492925](https://docs.msupply.org.nz/docs/cold_chain_equipment:configure?rev=1768492925)

Last update: **2026/01/15 16:02**

