

28.02. Remote Support - Installing TeamViewer

Over the years we have used a number of applications for providing remote support to clients. As of writing this (January 2020) TeamViewer is the software that we use.

By following the installation process below, mSupply Support will have **unattended access** to the computer. This means that mSupply Support will be able to access the server without user interaction.

If you do not wish for mSupply Support to have unattended access, you can install [TeamViewer here](#). If you are not using our customised installer, any support requests for mSupply to connect to a device will require a **TeamViewer ID** and **password**.

It is highly recommended that you **rename your computer** and **use our customised TeamViewer**, but it is not a requirement.

Installing the TeamViewer Host customised for mSupply

1. **Uninstall** - if TeamViewer is already installed on the machine, it is safest to uninstall it.
2. **Optional: Change the Computer Name** to something meaningful **before** installing TeamViewer. Using a name that reflects the site and / or user would be helpful. e.g. Auckland - Grafton Pharmacy - this is what the machine will show up as on our internal directory.




The TeamViewer Host module appears to only use the first **15** characters of the name. If you are connecting several computers, make sure that the first 15 characters are unique for each computer.

To change the computer name in Windows go to: **Control Panel > System and Security > System > Advanced system settings > Computer Name > Change...**




Computer description is NOT the Computer Name


3. **Download** and save the zip file at the following link
<https://www.msupply.org.nz/files/Teamviewer/TeamViewer%20Deployment%20December%202025.zip>
4. **Extract** the files from the TeamViewer Deployment .zip file onto the host computer.

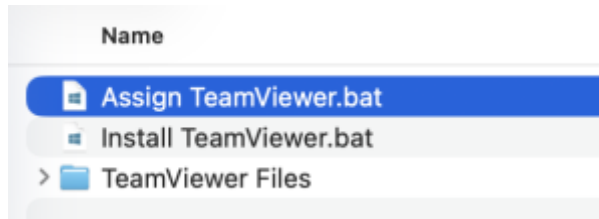

- 5. **Open** the folder you extracted
- 6. **Run as administrator - right click** on the file `Install TeamViewer.bat` and select **Run as administrator**. The script will check to see if it is running as an Administrator, if it isn't running as an Administrator, the installation will stop. If TeamViewer is already installed, you can try and use `Assign TeamViewer.bat` instead.

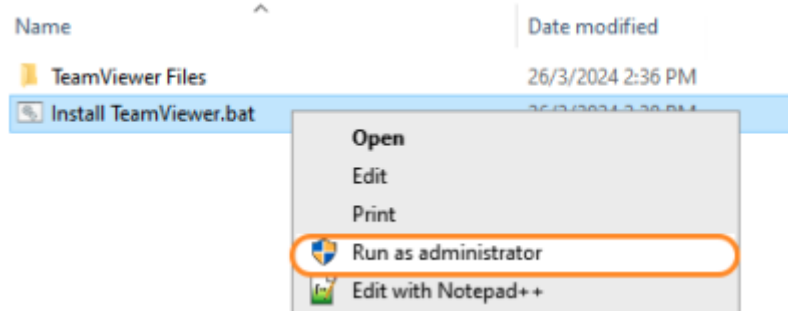

The Teamviewer Customised Host window should appear on your screen:

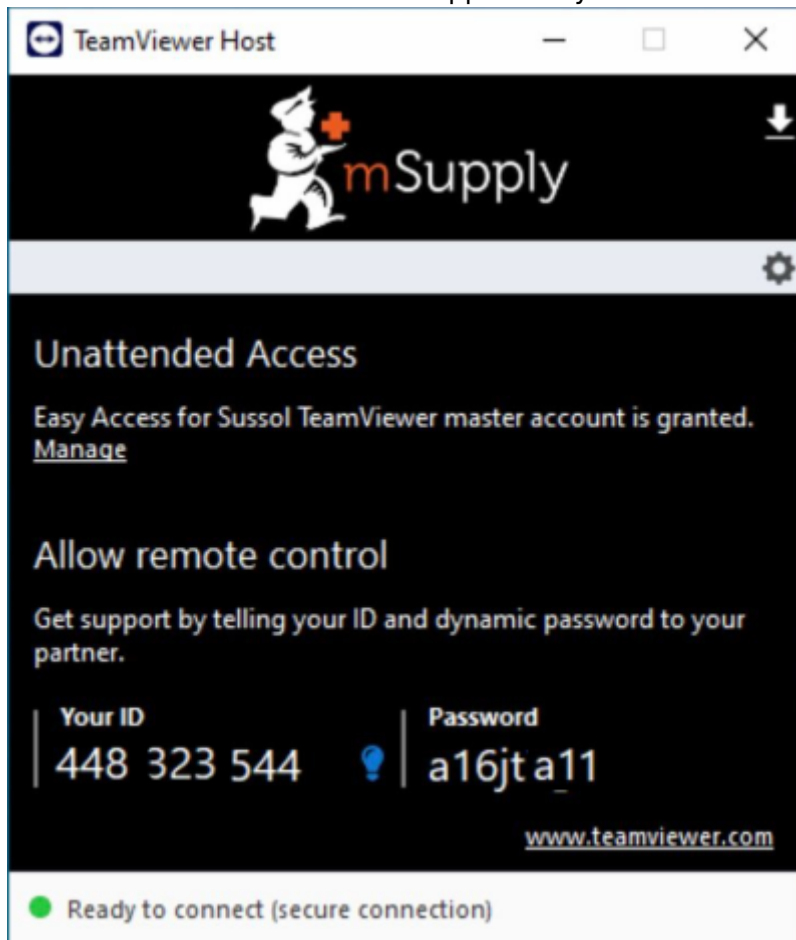

If the above window does not appear, the full installation has not completed. It definitely helps to run the batch (`Install TeamViewer.bat`) file with a reasonable internet



connection. If the internet connection is not good enough, then some or all of the actions in the batch file will not be complete. If that happens, please uninstall TeamViewer and re-run the batch file (as Admin) when the internet connection is better.

- 7. Notify mSupply Support** When you see this TeamViewer Host window, send us a message or email (on support@msupply.org.nz) and give us the **computer name** that you gave the machine in step 4 above. Otherwise, the TeamViewer ID should also work (a screenshot is fine). We will then be able to confirm for you that the remote access is operational and ready for us to provide you remote support.

From:
<https://docs.msupply.org.nz/docs/> - **mSupply documentation wiki**

Permanent link:
https://docs.msupply.org.nz/docs/faq:install_remote_support_software?rev=1775075764

Last update: **2026/04/01 20:36**

